


OPERATION AND INSTALLATION MANUAL

(MLG CF-680)

WIRELESS INTRAORAL CAMERA

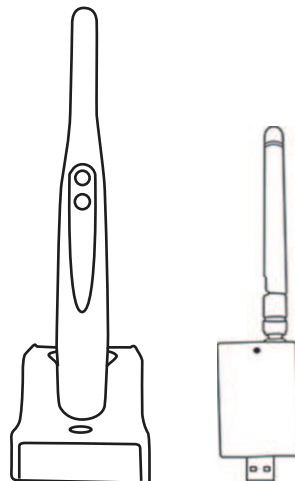

Product characteristics


1. Easy control: equipped with intelligent zoom button, free control, intelligent focusing, And no delay, zero interference, smooth operation experience.
2. Lens quality: multi-layer coating crystal lens, with dust and anti-fog function, can ensure the clear and stable shooting picture, even in complex oral environment can work better.
3. Flexible focus: Support 10mm close-range precision focus shooting, can also be 100mm Long-distance shooting to meet different observation needs.
4. Diverse connections: wireless charging is supported, convenient for daily use; It can be connected via USB The computer can also be wirelessly connected, and the data transmission method is flexible.
5. Clear imaging: equipped with CMOS chip, resolution of 13 million pixels, equipped with 8 white Color LED light source, Angle of view is 105°, output JPEG 1920 * 1080P format picture, Image quality is high.
6. Software adaptation: It can be used in combination with dental management software such as easy dentist, dental butler, and e-dental. AI computer software is also supported to facilitate data management and analysis.
7. System compatibility: can use Windows system (Win7/Win10)。

Product summary


parameter


Name: Intraoral camera	Voltage: 5V
Model: CF-680 II	Light: 8 pcs LED
Pixel: 1080P	Focus range: 10-100MM
Output: USB 2.0	Angle of view: 105°

PACKING LIST

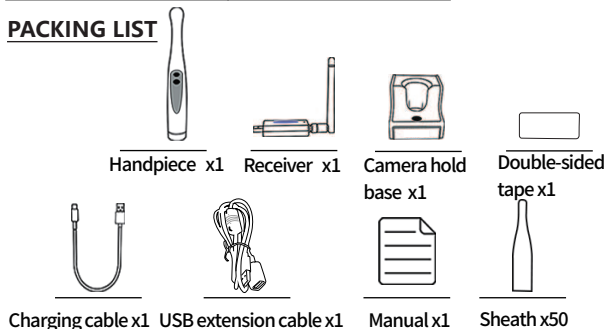

PRODUCT PROFILE

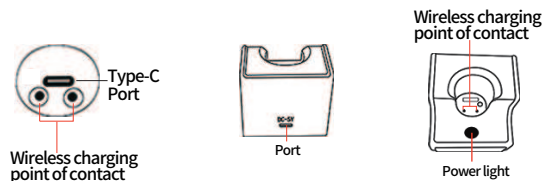

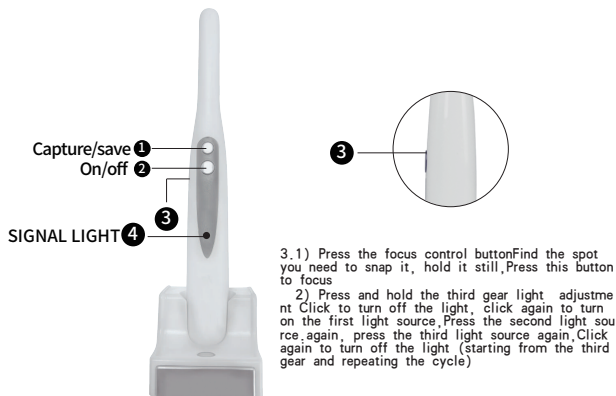

Operating


Computer conection




- 1.) Receiver connect to the computer
- 2.) Long press the power to turn on(notice: Press the power on button 2 seconds and the green light will lights up, blinking is already running)
(Automatic shutdown after 5 minutes of no operation)

Software installation


- 1: Connect the computer with the attached software security USB flash disk
- 2: Install the software to the computer
- 3: Double-click the software

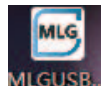

- 4: Click "User Management" to create a user profile

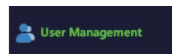


- 5: Click "Login" and then click "Video Capture Image Capture" to take photos


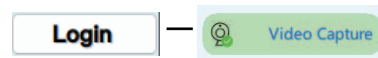


- 6: Enter the photo taking interface.
(All the latest pictures taken can be viewed at the top of the scroll bar)

- 1) Press the handle button to take photo
The latest picture will pop up at the top
- 2) Scroll up
Find the latest pictures taken

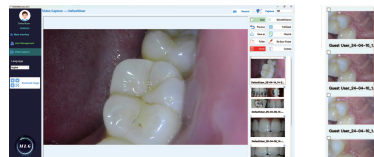

Toolbar description of the take photo:

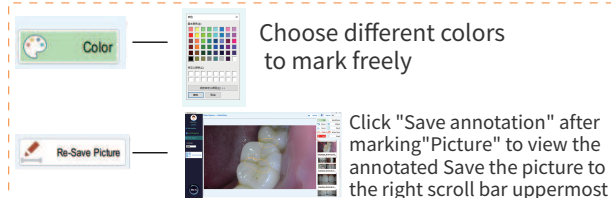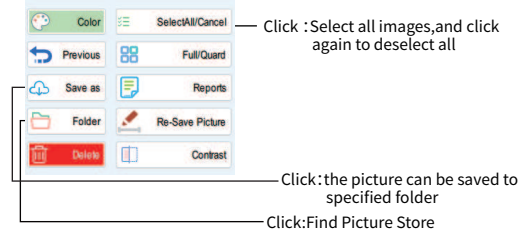




Back beat

Each one can be reshoot



Take a picture

Take four photos (take four photos to store)

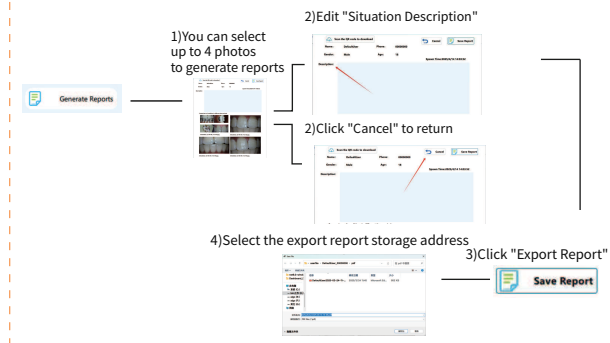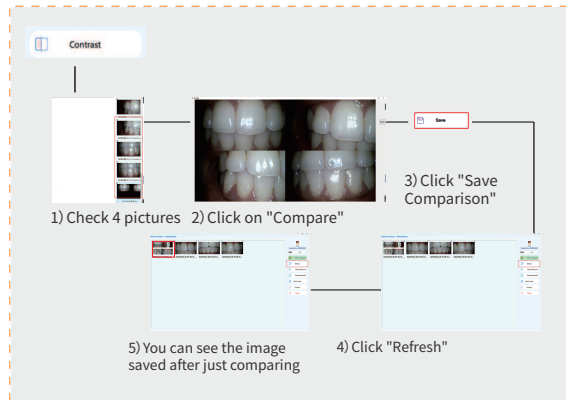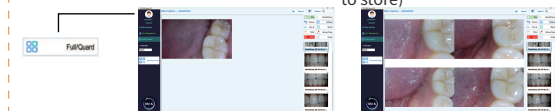

HOW TO SOLVE TROUBLES

- 1.the products can't be start-up when open (ON/OFF)switch) ?
 - 1).Check whether the battery of the camera has power, you can charge it first.
 - 2).After full charge still can not turn on, you can contact the dealer or manufacturer
- 2.computer can't reveal images?
 - 1) Check the connection between products and computer if it is correct
 - 2) Check whether the computer properties are set correctly
- 3.After turn on the lights are on but the connection is not successful?
When the signal light blinks, it indicates that the connection is successful. It takes about 3 seconds to wait for the CPU to run after the handle is started
- 4.in the process of use, the image will be a little blurry?
if you use it in the first time, it belongs to natural reaction. It will be better later on
- 5.the images of computer have spots?
 - 1) Check the screen if it has spot
 - 2) Check the lens if it has spot
 - 3) If the spot can't be cleaned, please contact the local dealer

NOTICE

Product

- Do not touch the signal connector with sharp objects,metal or liquid to avoid short circuit and product damage
- If you have not used this product for a period of time,please make it power off please fo not disassemble,which may damage the product

Power

- Please use the original manufacturer specified special power transformer, so as not to damage the product
- If you find the power cord or outlet is damaged ,please do not use it again

USB environment

- Please put the product in a dry and ventilated place
- Please do not put this product on a place are too cold, too hot or too humid , direct sunlight, dust and more places

Warranty service

Thank you for purchasing . In order to ensure our company's high-quality after-sales service, the following are instructions for warranty.

Product Warranty Conditions

- 1.The warranty period of the product is 1 year. During this period, if any problem of product, can inquiry the dealer or our company for consultation and maintenance.
- 2.The following cases are excluded from the warranty service, and we will charge a fee (charge standard = maintenance fee + replacement part fee + transportation fee or postage fee)
 - 1) The quality problems because of improper use or damage to follow product instructions
 - 2) When the product damage, replacement by consumers, dismantle it or any repair station without the company's approval.
 - 3) artificial damage, or damage caused by installation, movement, etc.
3. This product has passed the relevant quality management inspection as a qualified product before leaving the factory, but factors such as transportation cannot be ruled out. Once product breaks down, you can contact the dealer or our company with the warranty sheet to the warranty service

Contact alandental.com

Email: service.alandental@gmail.com

URL: <https://www.alandental.com>

Purchase date

Customer name

Telephone

Fault cause	Repair method	Date