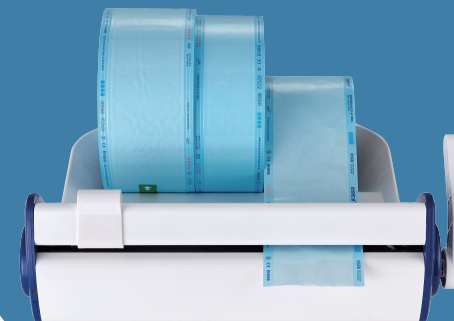




THE INSTRUCTIONS

- Please read the instruction manual carefully before using this product
- Please keep this instruction manual

⏪ Instructions

Please read the instructions carefully before using the sealing machine, safety tips and instructions, this manual will let you know the advantages about this machine, make you feel more comfortable when you using, please strictly abide by the operation, and reasonable maintenance, please keep instructions for reference.

If there are some problems when you using this machine, please contact the sales, we will provide the best service for you.

⏪ Avoid safety hazards

1. In the handling of sealing machine, should be handled gently, do not force too much.
2. When using the sealing machine, the sealing machine should always be facing up, and it is forbidden to put sideways or diagonally
3. Please do not splash water on the sealing machine in case of leakage
4. Please don't flip the machine

⏪ Product overview

1. Medical sterilizing bag is one of the sterilizing auxiliary products about sealing machine
2. It has the advantages of stable performance, beautiful appearance, convenient operation, control system, space saving, fault alarm, it will automatically power-off without operation for a long time, no assembly and so on. Built-in double-edge industrial cutting blade, sharp edge, flat cut, with strong wear resistance, and has long service life.
3. The power line is connected to 220V power supply and 60Hz power supply. Turn on the power switch, it takes about 3 minutes to heat the sealing machine, and the heating rod reaches 200±10°C. Automatic constant temperature, when the green light does not flash, it means the machine can be used normally.

⏪ Pay attention

The correct way to remove the sterilizing bag from the bottom of the bracket is nylon face up and paper face down (in contact with the heating bar)

⏪ Function is introduced

This product is an intelligent sealing machine. After electrification, the green light flashes and heats up until about 3 minutes later, when the green light stays on, it means the sterilizing bag can be carried out

Move the handle to the sealing (can be operated in both directions) pressure sealing paper, 1 of 3 seconds after the machine will drop with yellow lights flashing hint sealing effect has been achieved, can move cutter frame to carry on the prompt, if within half an hour without any operation, the machine will automatically power-off protection, need to hand back to normal work.

⏪ Technical parameters

Voltage	230V 60hz	indentation width	12mm
Rated power	100W	Weight	4.8kg
Sealing temperature	200°C±10°C	gross weight	7kg
Maximum sealing width	300mm	Package size	50*40*21cm

⏪ Equipment cleaning

It must be Power off before cleaning, it can use mild soapy water clean stains, do not let the liquid into the equipment, can use a soft cloth and a few of detergent wipe the plastic shell.

note

If there are any question, please contact the manufacturer or dealer immediately, do not privately disassemble and repair, it is man-made damage, the manufacturer is not responsible!