

Endo Smart +

New generation **Brushless Endo Motor**
Bring you a new experience in root canal treatment





6:1 Super mini contra-angle

Super mini contra-angle allows clearer operation field.


Combined length determination mode

Endo Smart + can be connected to the matched Apex locator to realize apex locating function.
By connecting the Endo Smart + to the matched Apex locator, the position of the file tip inside the canal can be monitored during the procedure and many automatic functions such as Auto Apical Reverse can be activated.

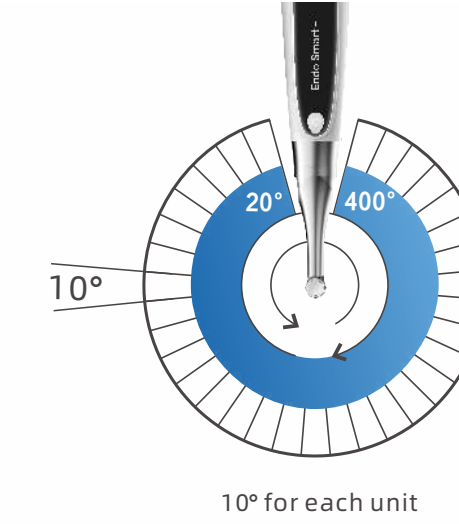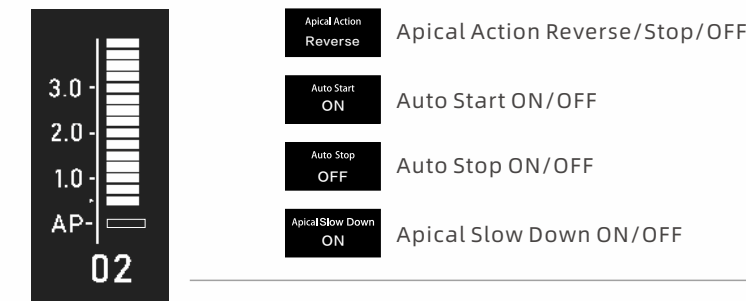

High performance brushless motor

- Adjustable reciprocating angle (10° interval; adjustment range: 20°- 400°); Accurate rotating angle reduces risk of instrument fracture and enables higher preparing efficiency.
- Super high torque, high speed, fast response, more accurate output, more stable, to ensure more efficient & safe preparation.
Speed : 100-1200RPM Torque : 0.4-5.0N.cm
- Smooth rotation, lower noise and small vibration in handpiece. Doctors feel good during operations. Patients feel more comfortable during operation.

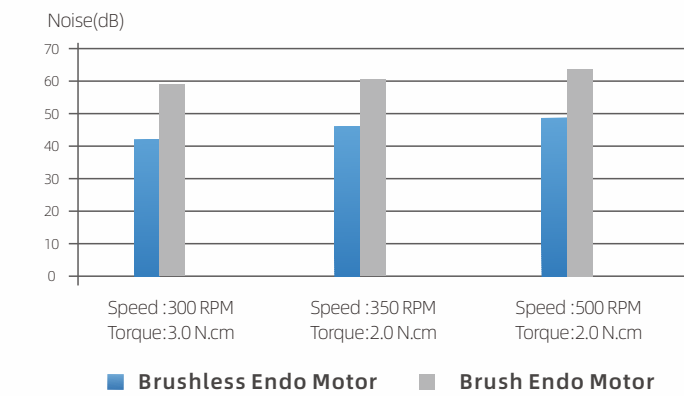

360° rotatable contra-angle



Three times longer service life



Compatible with all mainstream file systems



2000mAh Lithium battery



EXPORT PREMIUM QUALITY